




BIN CLEANING

M-LINE application

iwat[®]
a **TECNIPLAST** company
PHARMA

BIN CLEANING

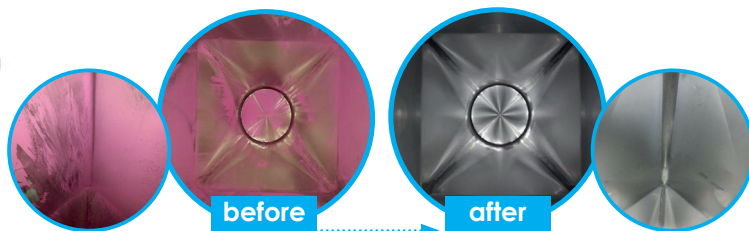

Do you have to...

- ▶ Automatically clean IBS of multiple sizes
- ▶ Minimise cross-contamination risks through an automated cleaning process
- ▶ Shrink storage space thanks to reduced downtime for cleaning
- ▶ Reduce water and detergent consumption

IWT have the perfect cleaning solution:

M-LINE HIGH PRESSURE WASHER

- ▶ One hydrokinetic lance to fit multiple size containers
- ▶ Lid to contain water splashes and allow cleaning of the manhole edge
- ▶ Non water-soluble sticky product: Results of cleaning validation
 - 2 min washing time with alkaline detergent at 0.1% concentration
 - 4 min rinsing with PW
 - 240 l / 63 gal water consumption
 - In-line conductivity control compliant
 - Lab TOC control compliant
- ▶ Lower water and detergent consumption: great savings in effluent treatments

to find out more information:

info@iwtpharma.com | www.iwtpharma.com