




**BIN
CLEANING**

Do you have to...

- Automatically clean IBS of multiple sizes
- Minimise cross-contamination risks through an automated cleaning process
- Shrink storage space thanks to reduced downtime for cleaning
 - Reduce water and detergent consumption

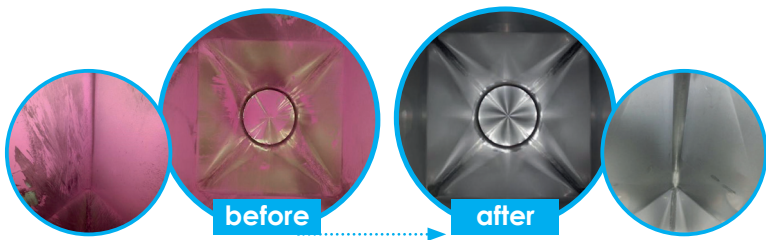

IWT have the perfect cleaning solution: **M-LINE HIGH PRESSURE WASHER**

- One hydrokinetic lance to fit multiple size containers
- Lid to contain water splashes and allow cleaning of the manhole edge
- Non water-soluble sticky product: Results of cleaning validation
 - ▶ 2 min washing time with alkaline detergent at 0.1% concentration
 - ▶ 4 min rinsing with PW
 - ▶ 240l/63 gal water consumption
 - ▶ In-line conductivity control compliant
 - ▶ Lab TOC control compliant
- Lower water and detergent consumption: great savings in effluent treatments

to find out more information: